




Anigen SYBR Green RT-PCR Kit

- *One tube use system for cDNA synthesis and PCR*
- *Applicable in commercially available Real-Time PCR apparatuses*
- *High Sensitivity and Specificity*
- *Perfect optimization for One-step RT-PCR*
- *Perfect RT-PCR possible by adding RNA template, primer, and RNase free water only*

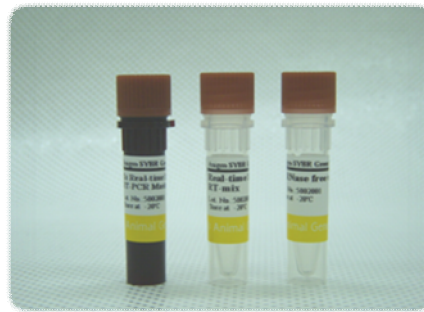

■ Description

Anigen SYBR Green RT-PCR Kit is workable for commercially available all real-time PCR apparatuses which detects SYBR Green I and it is available in most of RNA PCR samples.

It provides RTase, Hot start Taq DNA polymerase, SYBR Green, ROX reference dye, dNTP, and One-step RT-PCR buffer. It is only required that Template RNA and primer are mixed with the solutions of the Master Mix and RTase provided from the kit.

■ Specifications

- **Applications:** All Real-Time PCR
- **Cat. No:** RK50-02
- **Packing size:** for 1000 x 50 μ l Reaction
- **Contents**
 - 1 Bottle of 2 X concentrated of Real-time SYBR Green RT-PCR master mix(25m ℓ)
 - 1 Bottle of 0.5m ℓ of Real time SYBR Green RT Mix
 - 1 Bottle of 20m ℓ of RNase free water

